
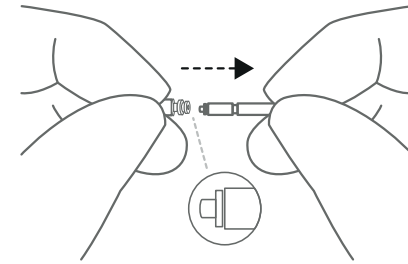
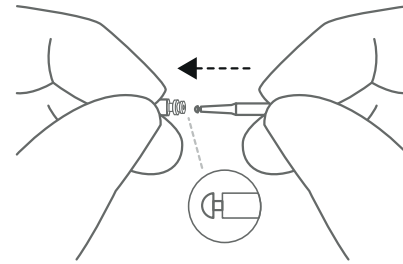



## How to change the wax guard

If you wear a dome, remove it before following this procedure. To replace wax guards:

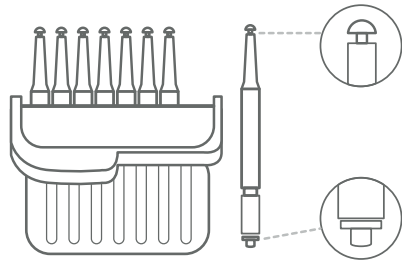


1. Carrying box with eight wax guard tools.
2. Insert the removal tip into the used wax guard until the shaft touches the rim of the wax guard.
3. Slowly pull the wax guard straight out.


The wax guard tool has two functions: a removal tip to collect the used filter, and a replacement tip with a white filter. To insert the new wax filters, follow these steps:

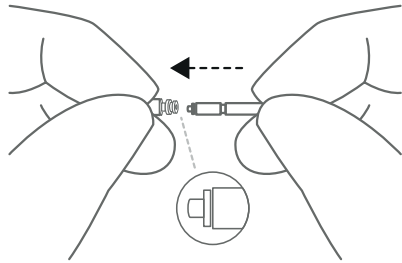

1. Insert the replacement tip of the tool into the sound outlet.

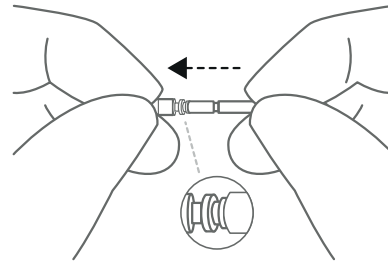

2. Gently press the replacement tip straight into the sound outlet until the outer ring lies flush with the sound outlet.

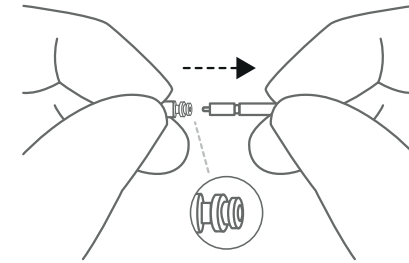

3. Pull the tool straight out - the new wax guard will remain in place. Please remember to re-attach the dome again, or attach a fresh dome.